

Dry Biochemical Analyzer

— Compass 2000



Low-maintenance



High precision



Disposable consumable



≤3min



4.3 Inch
touch screen



Constant temperature



Ultra portable



Optional printer

Test Parameters:

Lipids+Glucose

TC | TG | HDL-C | LDL-C | GLU

Kidney Function

Urea | Cre | UA

Liver Function

ALB | ALT | AST


Screening for Metabolic diseases

TC | UA | GLU


Screening of blood donors

Hb | ALT


Testing in three easy steps

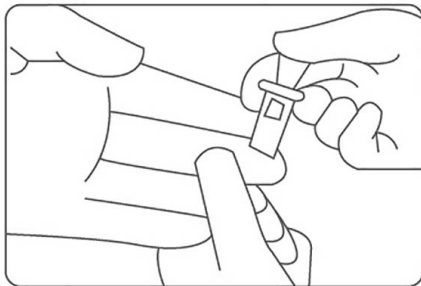

1. Puncture the finger using a lancet

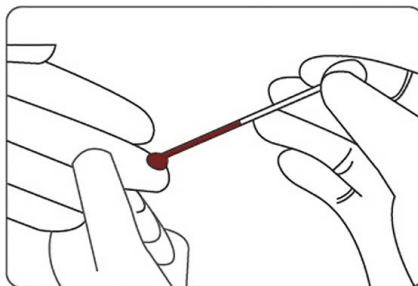

2. Draw a sufficient amount of blood when the blood drop is large enough.

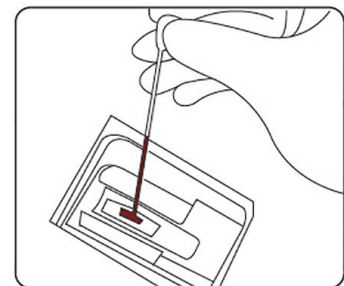

3. Drop the blood into the sample area of the test card for testing.

Specifications:

| | | | |
|------------------|--|-----------------------|------------------------------------|
| Detection method | Photometric detection | Measurement | 157.5*145.7*64.6mm |
| Detection time | ≤3min | Weight | 622g |
| CV | ≤5.0% | Data storage | 3000 results |
| Sample type | Serum, whole blood, plasma | Display | 4.3touch screen |
| Sample volume | 45μL | Storage environment | Temp: -20 °C ~ 60 °C; Rh:10% ~ 90% |
| Power supply | Adapter input AC100~240V,50/60Hz,30VA; output5V,2A Built-in Lithium battery: DC3.7±0.3V | Operation environment | Temp: 15°C ~ 35°C; Rh:≤85% |

Henan Forever Medical CO.,LTD

sales Office: Room 806, Floor 8 Wanda Center, Jinshui District, Zhengzhou, 450000, Henan, China

Factory: Industrial area, Changyuan Town, Xinxiang City, Henan Province, 450000, China

E-mail: forermed@yjcompany.com

Website: www.yongjiumed.com